

Revision: 201912
Runyes reserves the right to make modifications without prior notice.

DCI edge



DENTAL UNIT



DCI INTERNATIONAL
Add: 305N, SPRINGBROOK, RD. NEWBERG, OR 97132
Fax: 1-503-538-9302 Tel: 1-503-538-8343 <https://www.dciedge.com>

Runyes Medical Instrument Co., Ltd.

Add: No.456, Tonghui Road, Jiangbei Investment & Pioneering Park C, Ningbo
Tel: +86-574-27709922 Fax: +86-574-27709923
E-mail: runyes@runyes.com Website: en.runyes.com




Scan QR Code to
get informed on new product releases



A RELIABLE PARTNER

CENTRAL NEGATIVE-PRESSURE SYSTEMS

Dental chair with versatile and optimally designed features which provide great working ergonomics. Combination of doctor's convenient operation and patient's comforts, that is what we have been working on for many years. DCI dental unit, always at the forefront of the industry, keep researching and creating new developments, aiming to bring you a better control and sensational experience.

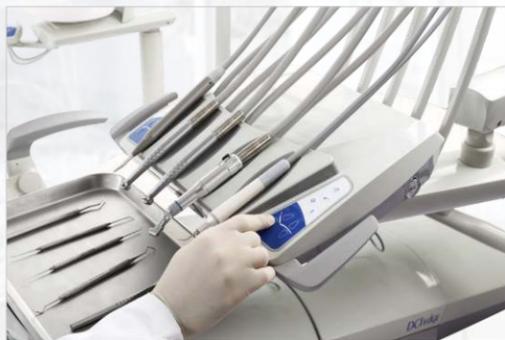

BEST IN QUALITY



DCI Edge

MAPLE 42U (UP MOUNTED)

The reliability of the DCI dental unit comes from high & accurate quality testing and solid materials. Original & imported components from DCI USA, international leading production facilities, and thousands of tests ensure the durability & fact of lasting a long time.

A combination of attractive appearance and dynamic design matches different styles of clinics. DCI dental unit meets the standards of ergonomics, which brings more comfortable experience to the doctors and the patients.

For different working mode, we optimize the operation process, providing you with the perfect solution. And our personalized settings give you the flexibility to configure your dental unit to improve efficiency.



This image is color:109



IT IS TIME TO CHOOSE A BETTER WORK PARTNER

You deserve the best, so do your treatment team and patients.



Chair Rotation
Up to 30°



Hydraulic Control
System



Dual Foot Panels



LED Operating
Light



Effective Disinfection
System



Flexible Choice
of Working Hand



Maple 42D(DOWN MOUNTED)

The cuspidor is made of anti-staining glass porcelain and can be removed for disinfection, with rotation of 90 degrees. It also has a dosing and rinsing function that can be programmed on the control panel.



Each pylon can be independently adjusted in angle, satisfying the habits of different doctors to achieve optimal efficiency.



Double pedal control, one for controlling the dental chair lift system and the other one for controlling the dental treatment system. A separate system is more stable and can avoid erroneous operations.



Prevent accidents by turning off the water and electricity main switch to eliminate safety hazards when not available at clinic.



Integrated catheter irrigation and disinfection system, can flush one by one or flush all tubings at the same time. Provides quick and effective cleaning and resist patient from cross contamination.



This image is color: I06

YOUR COMFORT IS OUR MISSION

DCI **EDGE**

MAPLE 62U
(UP MOUNTED)

Lowest height of 340mm from the ground, provides convenience to children and disabled patients while entry & exit.

The headrest can be adjusted up to nine angles, only need one hand to adjust the position rapidly and conveniently.




DCI edge

MAPLE 62U (UP MOUNTED)

USER-FRIENDLY TO MAKE WORK EASIER

Optimum work conditions for assistant makes your practice more productive. Since your assistant element is fully compatible with your dentist element, the chair, hygiene and spittoon related functions, are always easy to call up, using the direct selection keys.

Rotating instrument tray, can be switched around without any tools. So the dentists can share with the device whatever left-hand or right-hand.

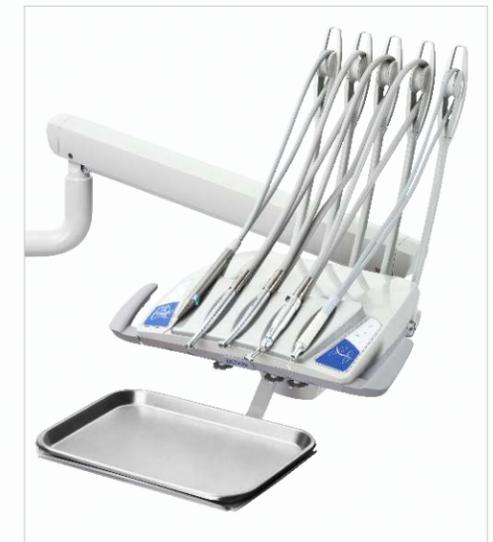

This image is color:I08

MORE ACCURATE MOTIONS WITH EASIER OPERATION



Ergonomic Design



Triple-Axis-Movable LED
Operating Light



Water Purification
System



Self-Cleaning Cuspidor




Safety Control Switch



Hydraulic Driver System



Balance arm

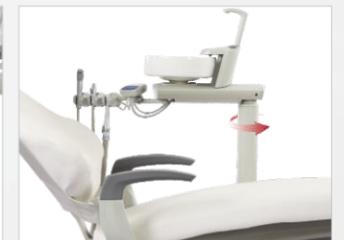

Assistant Set

Maple 62D(DOWN MOUNTED)



The ultra-thin upper narrow and wide backrest allows you to perform treatments in a comfortable position with your legs under the patient's body and the elbows to maintain their natural position, reducing fatigue and muscle tension.

Patented electro-hydraulic drive system, with smooth start and static setting, load weight of more than 220 kg, lifting speed can be manually adjusted.

This image is color:I03



Triaxle positioning adjustment lamp head, flexible positioning, always aiming at the patient's mouth, providing you with a high quality working light source.

The assistant position rotates around the dental chair, allowing for both left-handed and right-handed operation.

The left and right air-balanced arm passes under the dental chair and does not cause the balance arm to pass over the patient and to make the patient feel tense.



LED oral light, which can be inductive control or manual control also have the function of blue light cut off. Available in white and yellow light sources, which can protect the curing material from abnormal curing.



The rotatable and adjustable left and right armrests can be rotate to a position which will not interfere with operation, so it is much easier for the dentist to operate and easier for the patients to get up and lay down.

Equipped with locking lever so that dental unit seat can be rotated 30° to the left or right, Which is more convenient to cooperate with the surrounding medical equipment.





Built-in water bottle



Power box



Hidden host



Electronic clock

DCI edge
Maple 82



- Side cabinet design is brief but not simple.
- Self-contained water supply system, the instrument tray can be installed inside the side cabinet

*DCI*EDGE

