

Woodpecker Brushless Endo Motor

Endo Pace

MEDIFIX
ELECTROMEDICAL SYSTEMS



Endo Pace

Woodpecker Brushless Endo Motor

- **Lightweight handpiece**
Easy preparation
- **Brushless motor**
Strong power
- **10 times the service life of brush motor**




High performance brushless motor

High torque and fast response

More accurate and more stable output,
enabling more efficient preparation

Speed: 100rpm-1500rpm


Torque: 0.4N.cm-4.2N.cm

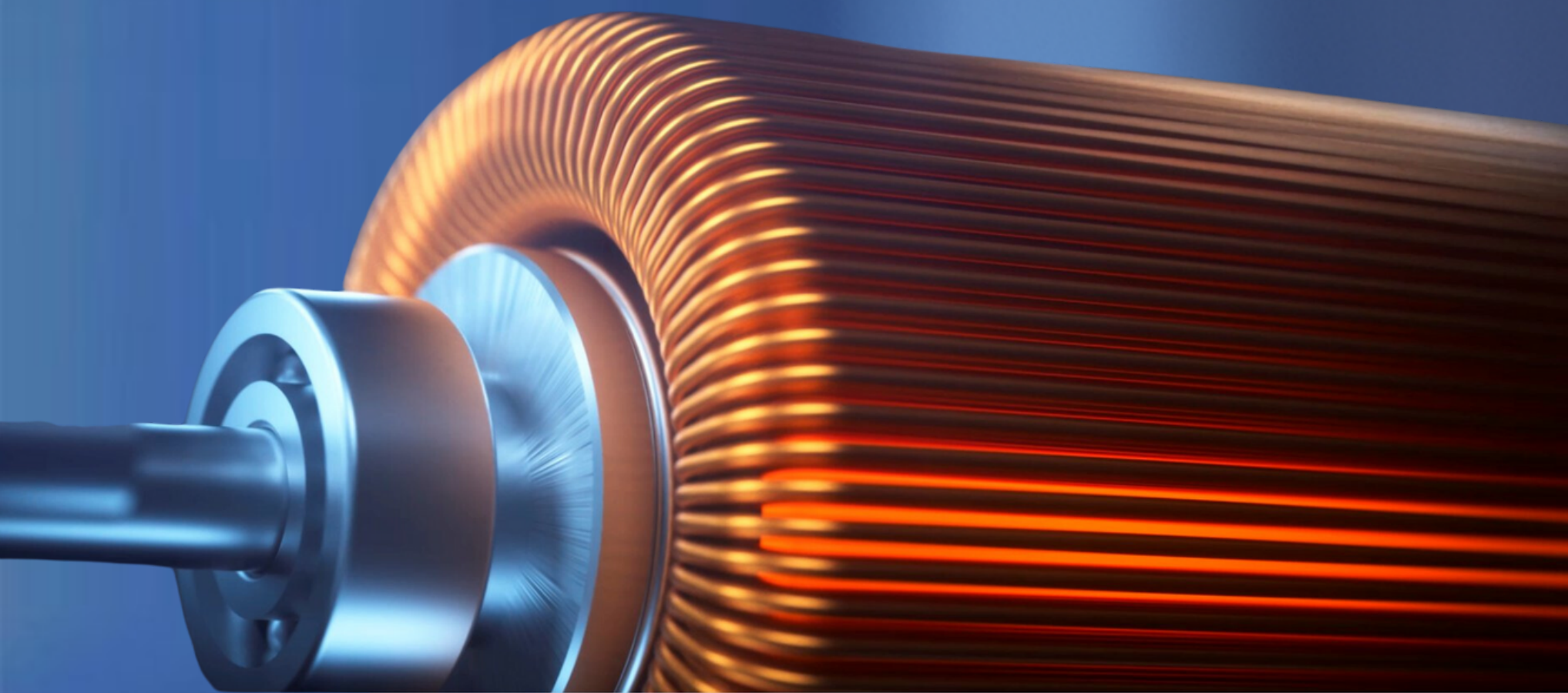

High performance brushless motor for low wear, low heat and longer service life

10 times the service life of brush motor.

Smoother rotation and lower noise.
Stable cutting can still be ensured at
the narrow or curved apex.

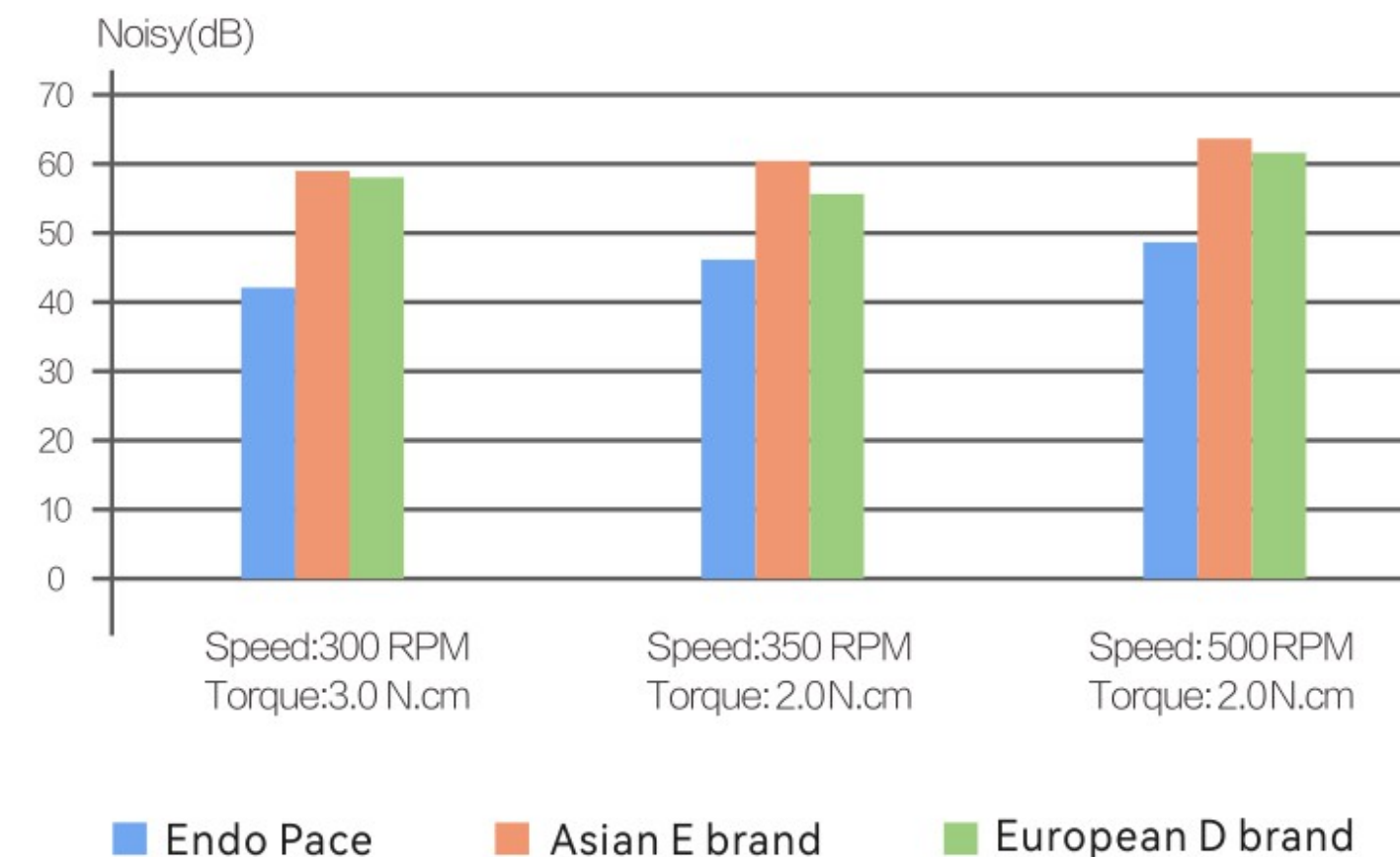


Adjustable reciprocating angle with 10° interval

Adjustment range of forward and reverse

reciprocating angle: 20°- 340°

Accurate control of the rotation angle reduces the risk of instrument separation and enables more efficient and safer preparation.

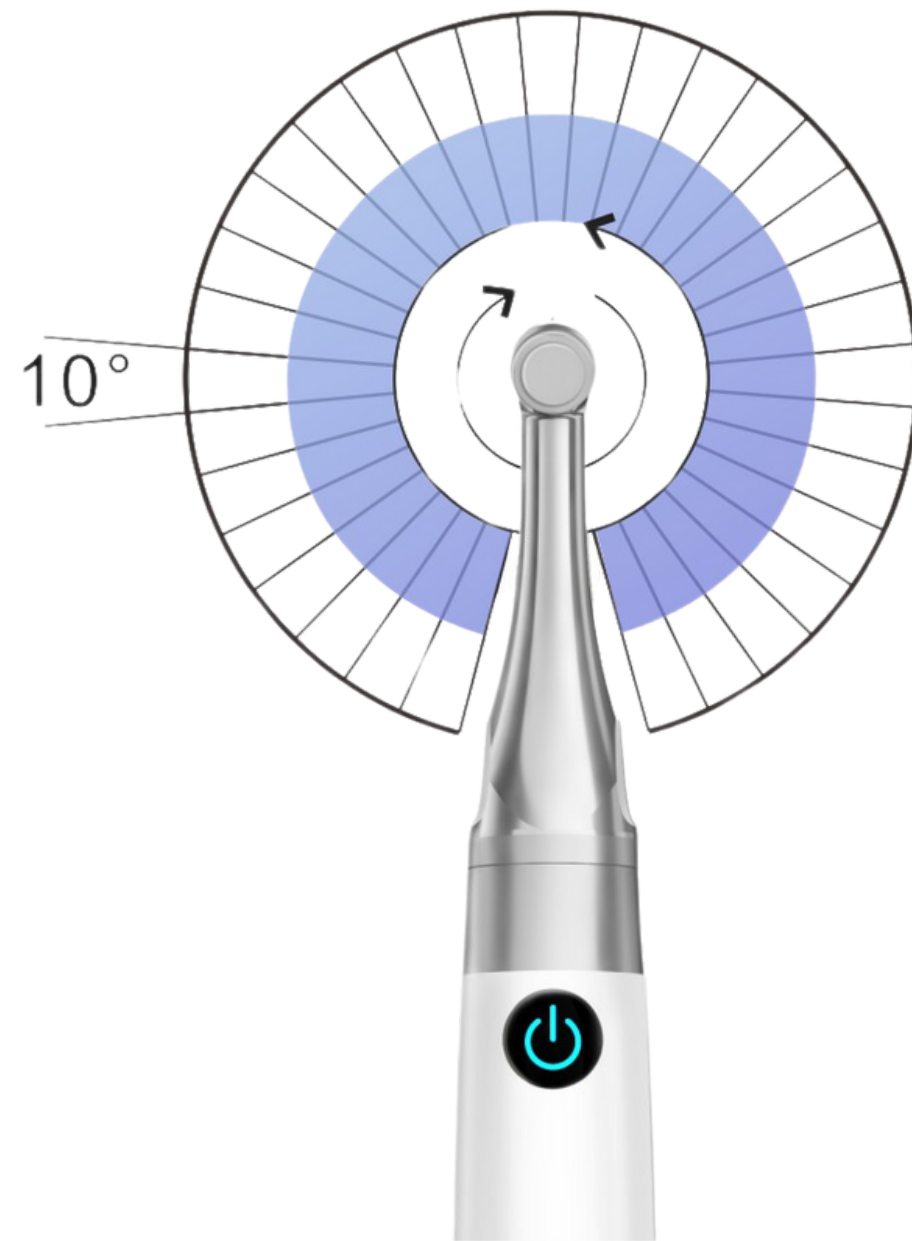

Lightweight cordless handpiece offers better hand feeling

Weight: 106.8g (handpiece+contra-angle)


Control the handpiece with wireless signal, without the restraint of tail cord

Well-balanced design allows easier preparation



360°rotatable contra-angle

More flexible operation,
facilitating the preparation of teeth
at diferent positions

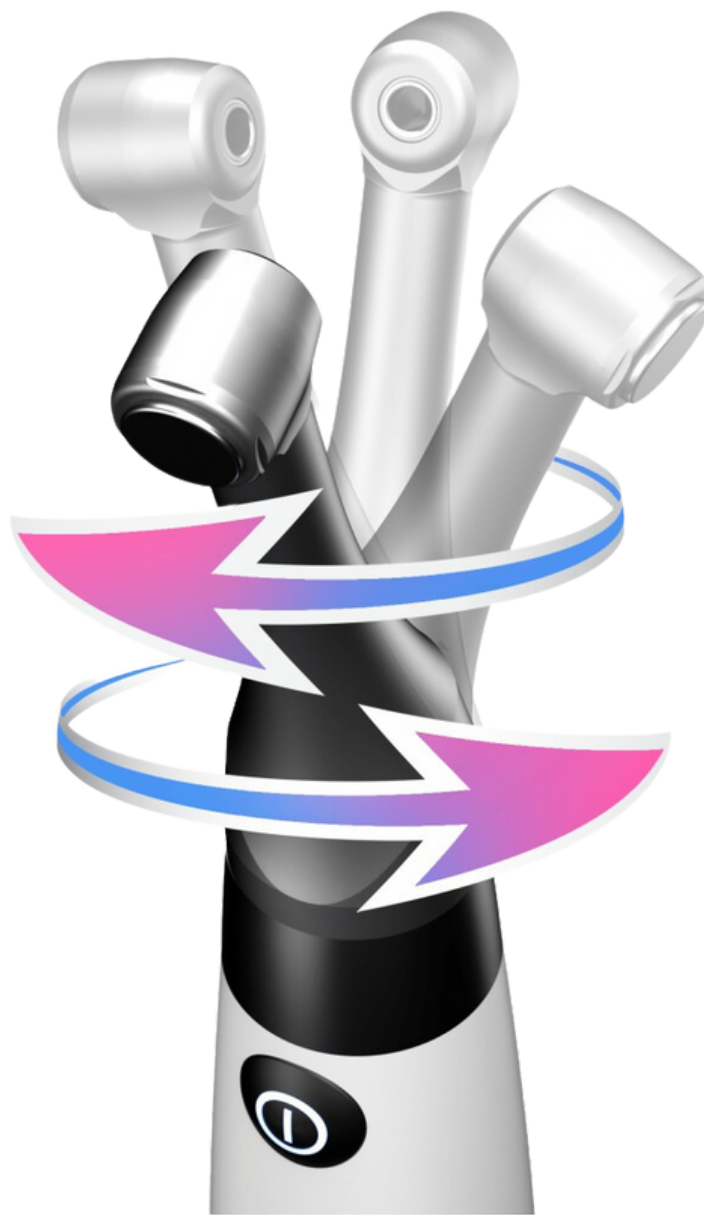

6:1 mini contra-angle delivers broader operation field

Compact but powerful
6:1 speed ratio ensures low speed
and stable rotation under high torque

Head height: 9.7mm
Head diameter: 8mm

Neck diameter: 6mm




**Can be connected to
an Apex Locator**



**4 working modes
guarantees safer and
more efficient preparation**

Forward rotation mode / Reverse rotation mode
Reciprocating mode / ATR mode

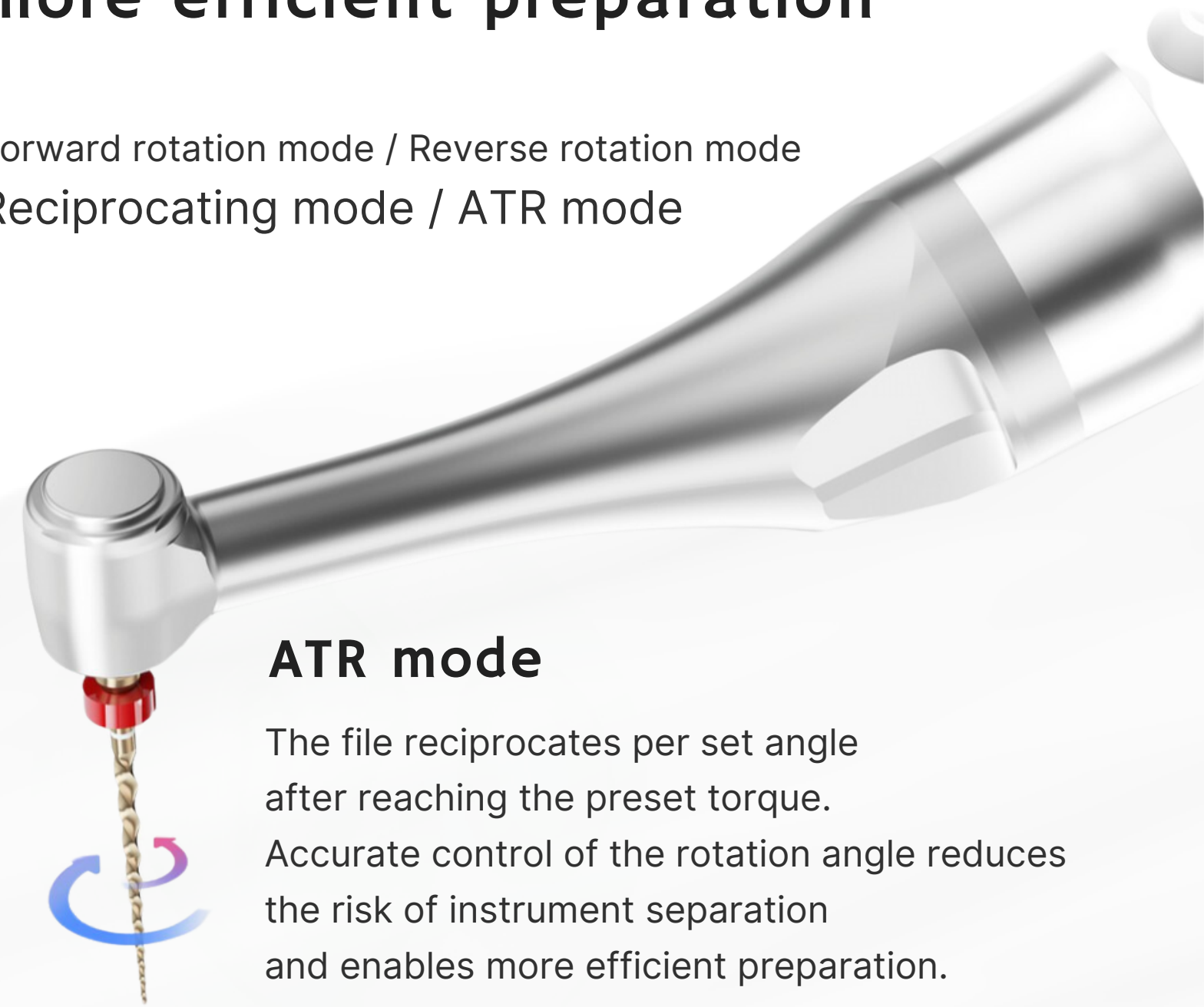

ATR mode

The file reciprocates per set angle after reaching the preset torque. Accurate control of the rotation angle reduces the risk of instrument separation and enables more efficient preparation.



 8626 4641

 8335 8987

 ventas@medifixcr.com