




# WOODPECKER HP SERIES DENTAL TURBINE HANDPIECE

Rebirth in Dentistry  
Making Smiles Shine Brighter



Model	HP33		HP22		HP11		HP13	
	NL, KL	M4, B2, N, K	NL, KL	M4, B2, N, K	NL, KL	M4, B2, N, K	NL, KL	M4, B2
Light	Optical Fiber	/	Optical Fiber	/	Optical Fiber	/	Optical Fiber	/
Chuck	Push-button		Push-button		Push-button		Push-button	
Chuck	Φ12 x13.1 mm		Φ11.2 x13.1 mm		Φ10.2 x12.1 mm		Φ12 x13.2mm	
Air Pressure	2.5-3.0kgf/cm <sup>2</sup>		2.2-2.7kgf/cm <sup>2</sup>		2.0-2.5kgf/cm <sup>2</sup>		2.5-3.0kgf/cm <sup>2</sup>	
Water Pressure	0.08-0.2MPa		0.08-0.2MPa		0.08-0.2MPa		0.08-0.2MPa	
Speed	280000-450000r/min		280000-450000r/min		280000-450000r/min		280000-450000r/min	


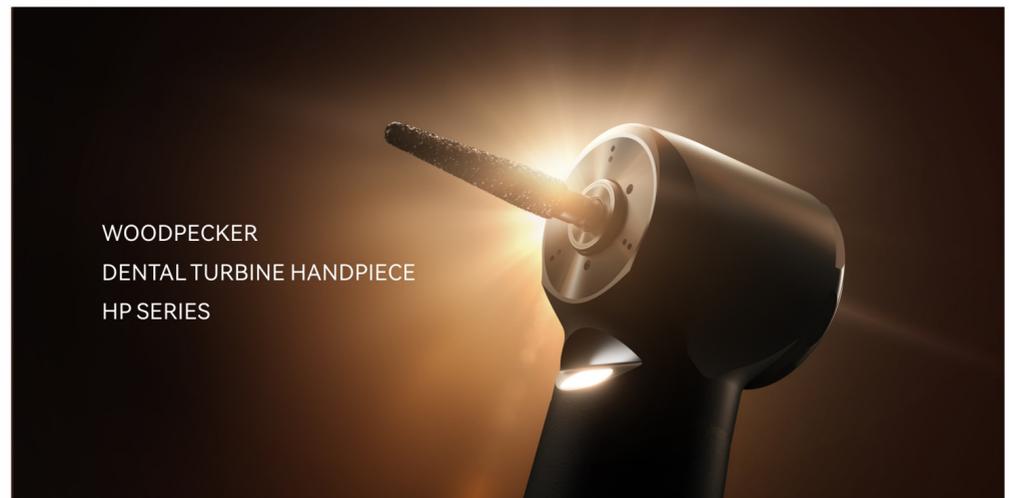


WOODPECKER  
DENTAL TURBINE HANDPIECE  
HP SERIES




Guilin Woodpecker Medical Instrument Co.,Ltd.  
Information Industrial Park, Guilin National High-Tech  
Zone, Guilin, Guangxi, 541004 P.R. China  
Sales Dept.: +86-773-5873196/2350599  
After-sales Service Dept.: +86-0773-5827898  
E-mail: woodpecker4@glwoodpecker.com  
Website: http://www.glwoodpecker.com

DISTRIBUTED BY:

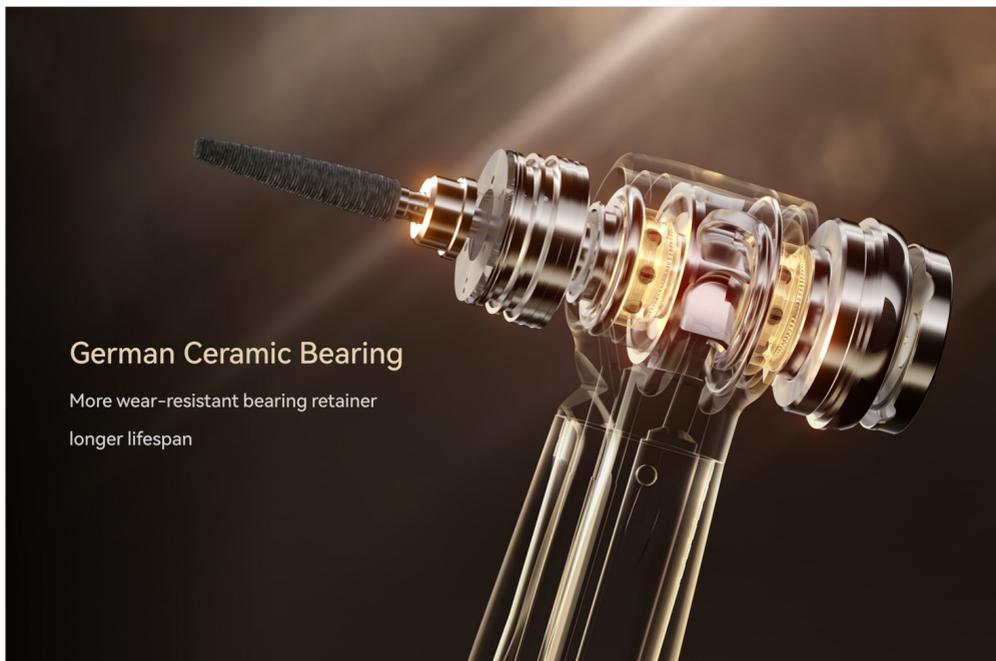

## German Ceramic Bearing

More wear-resistant bearing retainer  
longer lifespan

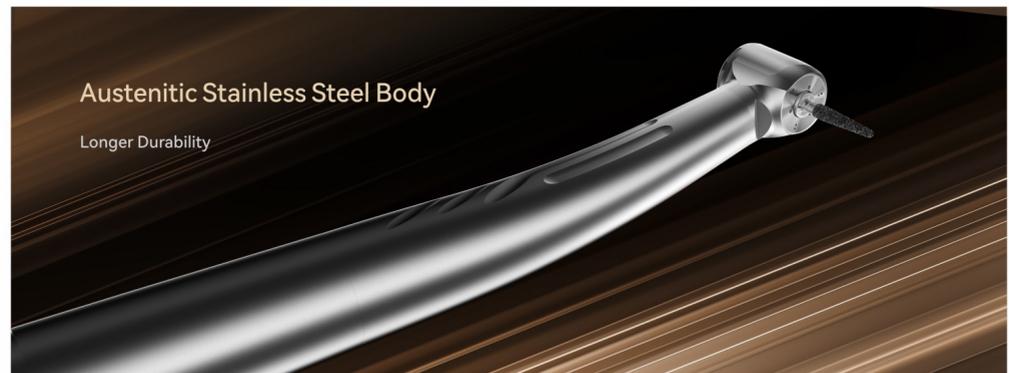


## Austenitic Stainless Steel Body

Longer Durability



## Well Sealed Clamping Structure

World's first patented technology  
avoids cross-infection

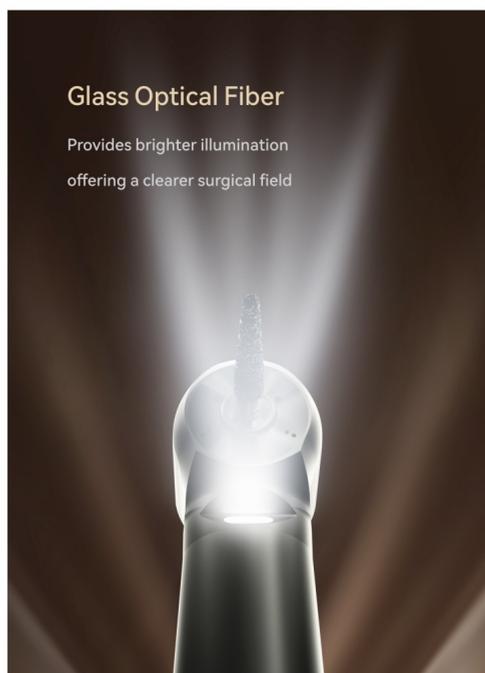

## Glass Optical Fiber

Provides brighter illumination  
offering a clearer surgical field



## Triple Spray Water-air separation

360° all-around  
cooling



## Striped Pattern

Comfortable grip  
easy to clean



## Thickened front end edge of the head

Better impact resistance  
reducing the risk of damage