

Woodpecker Ultrasonic Scaler

UDS-A LED



MEDIFIX
ELECTROMEDICAL SYSTEMS

A good helper for dentist

Optical handpiece, more convenient for clinical operation.

Automatic frequency tracking ensures that the machine always works on the best frequency and performs more steadily.

The handpiece is detachable and can be autoclaved to the high temperature of 135 and the pressure of 0.22MPa.

Digitally controlled, easy operation and more efficient for scaling.




Detachable handpiece with light

Shell of the handpiece is made of special imported plastic from America.

Without cracking or deformation under the pressure.

Detachable handpiece with light, better visibility.

Note: This type can match with handpiece without light.

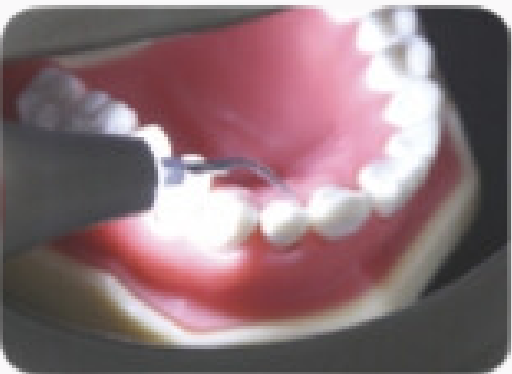

LED with light effect




Without light efect



NITI-File(included)

15# 20# 25# 30#

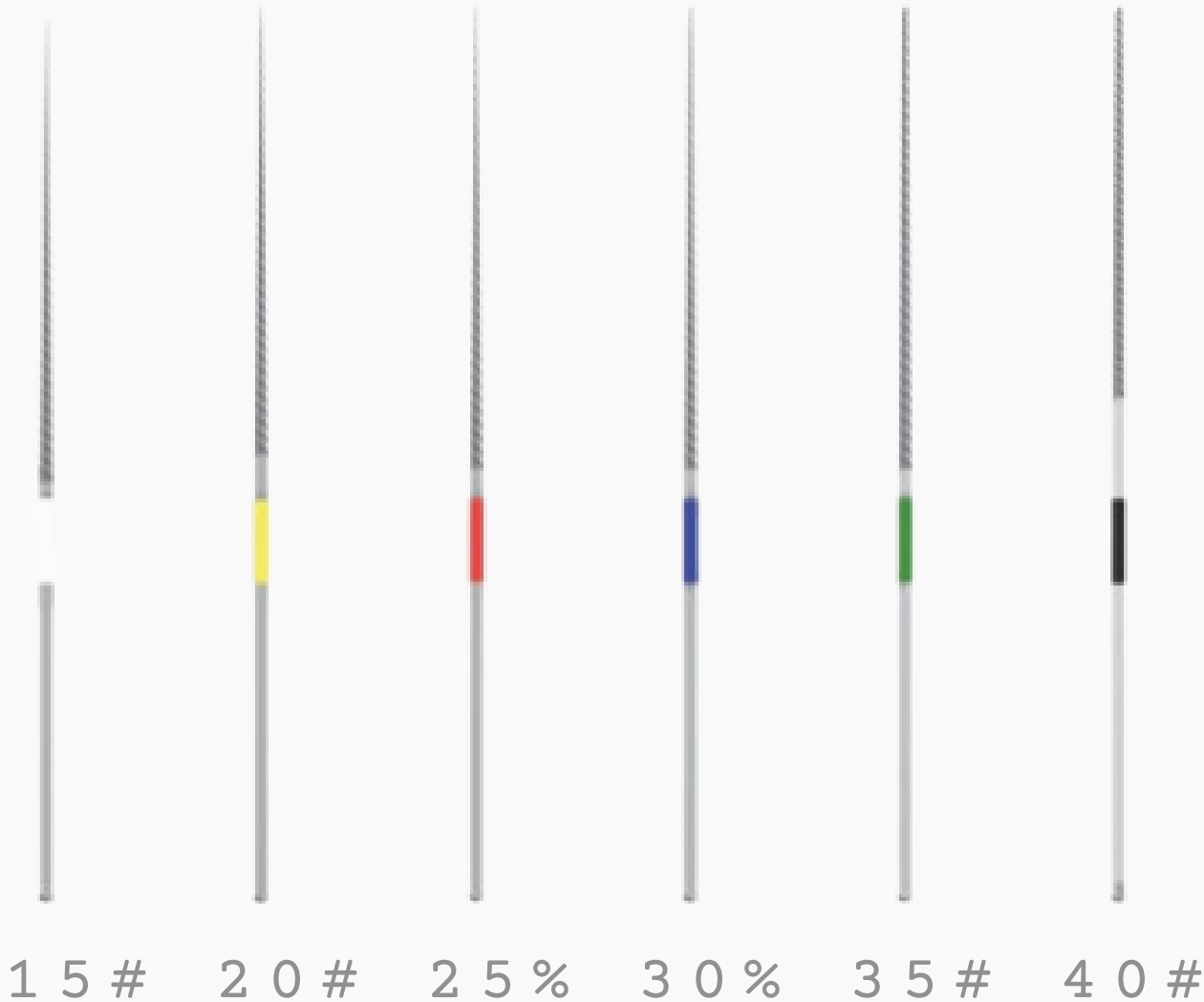

TIPS

> Supragingival scaling



> Implant maintenance



P90

Implant maintenance
Clean implants and restorations

P90

Subgingival scaling



Root canal treatment

Endo scaling



E1

E2

For removing the calcification and bad fillings in root canal



E3D

For removing the calcification in the coronal third of the root canal



P4D

Root canal retreatment, removing foreign body



(E6 E7)

Special material (E6, E7) , extremely thin and soft, good conductivity to ultrasound



E6

For removing foreign body in 1/3 deep root canal



E7

For removing foreign body in 1/3 root canal



E14

For removing foreign body and broken instruments in root canal (stainless steel)



E15

For removing root canal fillings and broken instruments (stainless steel)


Looking for root canal orifice



E12D

Diamond-coated, ball head, for looking for root canal orifice with security, accuracy and efficiency

Cavity preparation

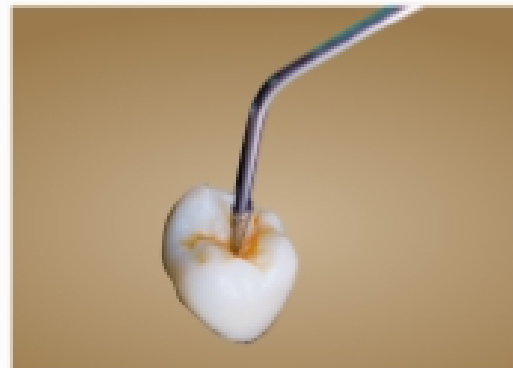


G30

For the preparation between pits and fissures



G31

For the preparation of mesioocclusion cavity

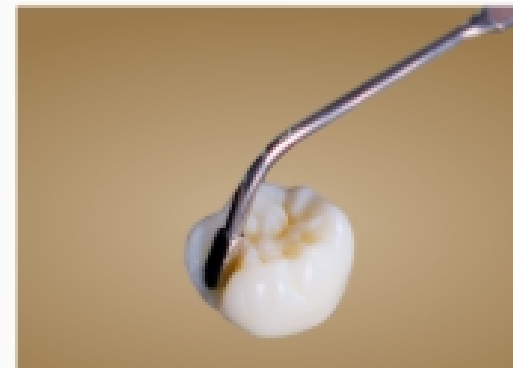

G32

For the preparation of distal mesial cavity



G30



G31



G32



SB1



SB2



SB3



 8626 4641

 8335 8987

 ventas@medifixcr.com